



BabyCare

Pocket doppler

BabyCare is a compact pocket doppler, which provides crystal clear fetal heartbeat sound. Being operated with battery, it allows babies to be checked anywhere, at any time. Using high sensitive and harmless ultrasound sensor, it detects and displays the accurate FHR on the LCD. BabyCare is especially suitable for routine use in OB/GYN. CE, CSA approved..

Features

- Accurate FHR detection & display
- Clear sound offered from separate speaker
- Cap-shaped speaker for probe protection
- Compact size and battery operated
- Harmless ultrasound sensor
- Easy-to-use







Maximized Convenience for FHR detection in OB/GYN

Technical Specification

| | |
|------------------------|-------------------------|
| Display | 3-digit segment |
| FHR Measuring Range | 50-240 bpm |
| Measurement Resolution | ± 3 bpm |
| Speaker output power | 1 W |
| Ultrasound Frequency | 2.0 Mhz |
| Power | 9V battery |
| Sensitivity | 12 Weeks Fetus |
| Weight | 290 g (without battery) |
| Standard Accessory | Gel, 9V battery |

This specification is subject to change without notice.

Quick operation

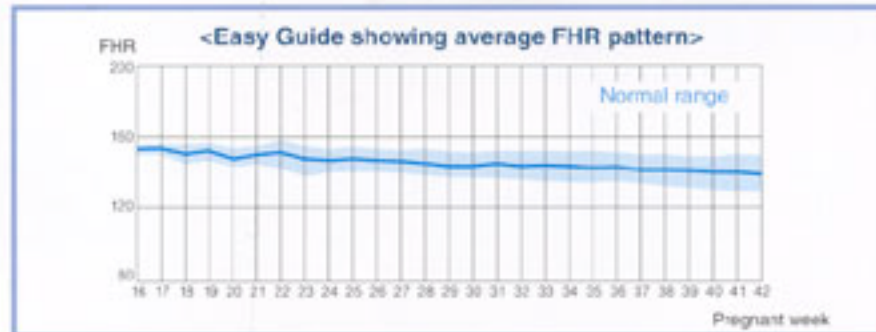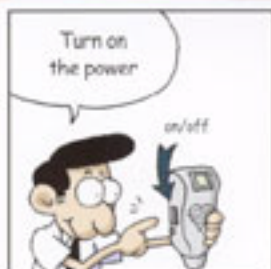

3F DAEYANG B/D., 999, DAECHI-DONG,
KANGNAM-GU, SEOUL, KOREA

Tel : +82 2 3468 3642 / +82 2 3468 3650

Fax : +82 2 3468 3699

E-mail: cichoi@bio2net.com / hwack@bio2net.com

www.bio2net.com www.agamom.com